# Artist of Science

## **DW-LDO-101 Series Forced Air Drying Oven**



- Horizontal forced air type good temperature uniformity
- Microprocessor controller (with temperature correction and timing function), accurate and reliable.
- Large LCD screen display.
- High quality stainless steel chamber, removable shelf, easy-to-clean.
- Equipped with multiple sets of heaters.(For DW-LDO-101 Series)
- Adopt temperature selection switch which can select the grade according to the heating speed and operating temperature.
- (For DW-LDO-101 Series)
- Silicon sealing ring for reliable sealing.
- Equipped with leakage protection.
- Equipped with spare temperature control. Ensuries to work normally when the main temp.control failed.
- Anti-hot handle.


#### **Technical Parameters**

- Temp. Fluctuation:±1°C.
- Temp. Uniformity: $\pm 2.5^{\circ}C\%$ .
- Temp. Resolution: 0.1°C

#### Optional Parts

- Multi-segment programmable control
- Built-in Printer
- RS485 interface



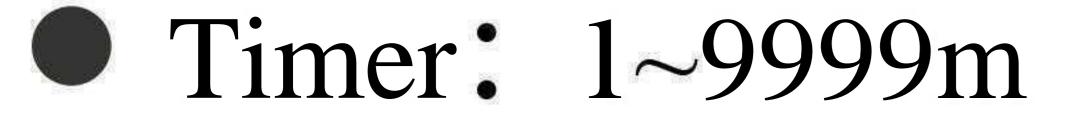

### Specifications

Model	Volume(L)	Chamber Size (WxDxH)cm	Power Rating (W)	Temp. Range (°C)	Shelf/Max	Net/Gross(kg)	<b>Power Supply</b>
DW-LDO-IOI-0	43	35×35×35	<b>1300</b>	RT+5~300	2/4	40/60	AC 220V 50HZ
<b>DW-LDO-101-1</b>	81	45x40×45	2000		2/6	55/85	
<b>DW-LDO-101-2</b>	138	50×50×55	<b>2500</b>		3/7	<b>80/115</b>	
<b>DW-LDO-101-3</b>	252	60×60×70	<b>3500</b>		3/10	110/150	
<b>DW-LDO-101-4</b>	<b>640</b>	80×80×100	5500		4/15	<b>170/250</b>	AC 380V 50HZ
<b>DW-LDO-101-5</b>	<b>960</b>	100×80×120	<b>7000</b>		5/18	250/340	

