

Horizontal Laminar Flow Cabinet

Model SW-CJ-1CU

Working Principle

First, the air is inhaled through the fan and then filtered by the HEPA filter through the static pressure box. After that the filtered clean air is sent out in the state of vertical or horizontal air flow, so that the operation area continues to reach the 100-level cleanliness under the control of clean air and to ensure the production of environmental cleanliness requirements.


Application Range

Horizontal flow work bench noise is compared small, so it is usually used in electronic industry, which has little effect to health. Horizontal laminar flow is widely used in electronics, national defense, precision instruments, meters, pharmaceutical industries, etc.

Features

- 1. Adopting arbitrary positioning sliding door system
- 2. The shell is welded by cold-rolled plates, the surface is electrostatically sprayed, and the work surface is SUS304 brushed stainless steel, which is corrosion-resistant and easy to clean.
- 3. Digital display LCD control interface, more user-friendly design
- 4. Equipped with HEPA air filter, with primary filter for preliminary filtration, which can effectively extend the HEPA filter life span



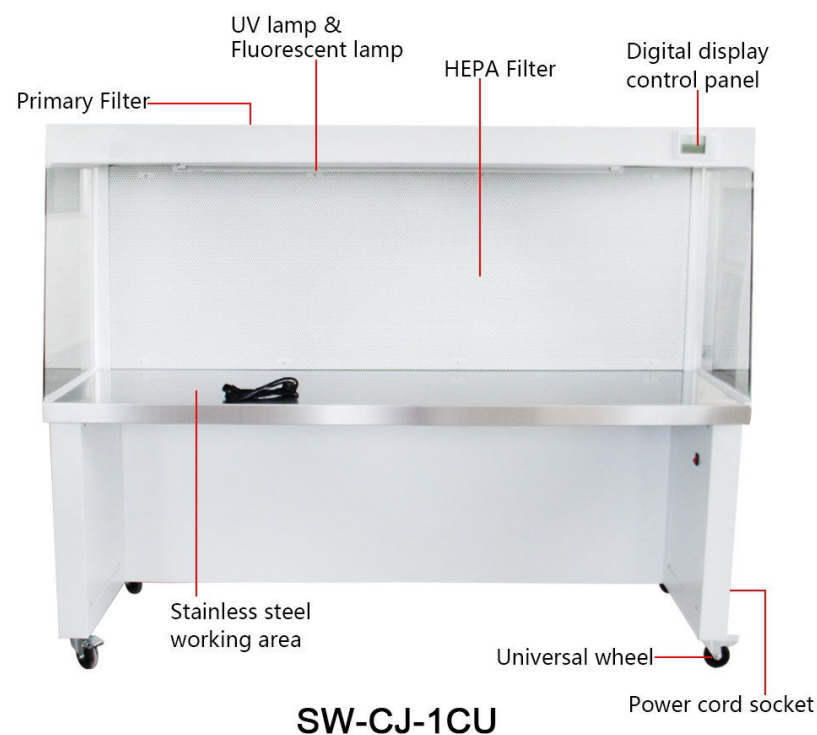

SW-CJ-1CU

Parameters

Model	SW-CJ-1CU (opening type)
Clean grade	Grade 100(209 E U.S. Federal)
Number of bacteria	≤0.5per utensil.hour(90mm utensil)
Noise	≤62dB
Power	220V/50HZ
Wind speed	0.25-0.45m/s (slow,medium and fast 3 speeds)
illumination	≥300LX
Working area dimension(W*D*H)	1770x510x600mm
Overall dimension(W*D*H)	1800x740x1440mm
HEPA filter size and quantity	1700x560x50x 2
Fluorescent/UV lamp specifications and quantity	20Wx 36Wx 2
Operator	Two persons