


# MPC-5V1500

## 2~8 °C Pharmacy Refrigerator



**MPC-5V1500** 2~8°C Pharmacy Refrigerator (1500 liters) is widely used in hospitals pharmacy stores and laboratories. This is for placing and preserving storage such as drugs, vaccines, enzymes, hormones, stem cells, extracted RNA and gene libraries and important biological and chemical reagents.

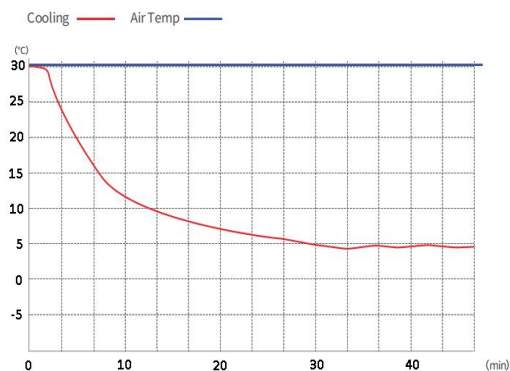

MPC-5V1500 PULL DOWN TEST

### FEATURERS



#### Fast pull down

Real air-forced cooling system



#### Super large space

1500-liter (53 cu.ft) space, More storage



#### Suction circulation system

Big fans with less temp. fluctuation & better



#### Universal castor wheels

5 proven heavy-duty castors



#### Lockable doors

3 Standard locks with keys for 3 doors



#### SUS 304 stainless steel

Internal & external both high quality SUS304



#### Precise temp. control

Temp. control accurate to 0.1°C



#### 5 types of alarm

Audible and visual alarms



#### Variable space

Height adjustable shelves



#### Eco-friendly

CFC-free refrigerant and energy-saving



## STRUCTURE

High-precision  
control system



Lock&keys





3 Large fans  
Internal LED tube



Optional Accessories

Data logger  
USB interface  
Chart recorder

Lockable  
castors

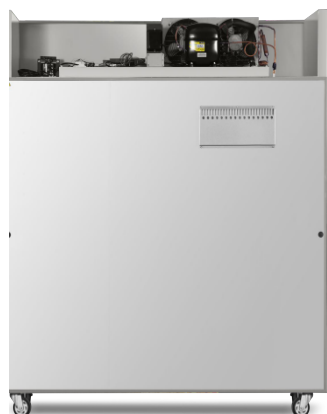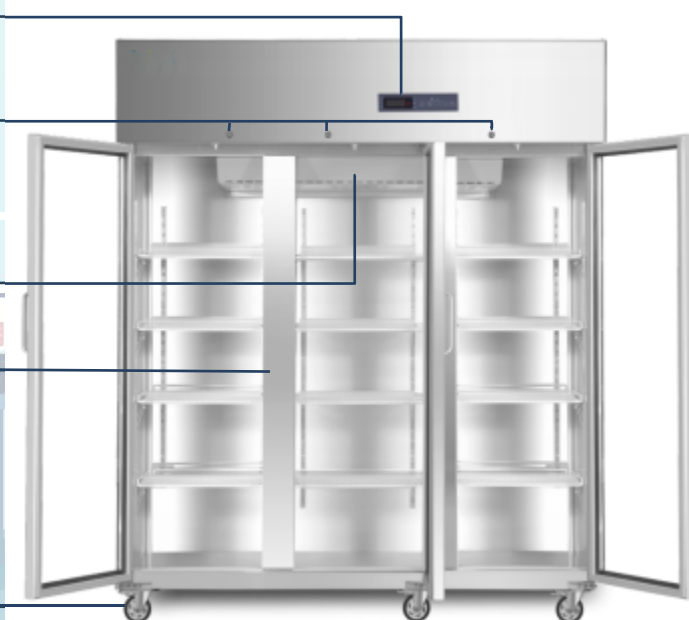

Imported brand compressor

High/low temp. alarm

Power failure alarm

Thermostat failure alarm

Door ajar alarm



Real air-force circulation and excellent refrigeration system pull down the chamber temperature from 25°C to 5°C in 35 minutes (at ambient temp. 25°C). And keep the chamber with good temperature uniformity and fluctuation.

## PARAMETER

Performance		Refrigeration System		Alarms	
Cooling method	Forced Air Cooling	Refrigerant	R134a	High/low temperature	Y
Defrost	Auto	Temp. controller	Microprocessor	Power failure	Y
Shelf/drawer	Shelves/12	Sensor	NTC	Sensor error	Y
Temperature range	2°C~8°C	Display	Digital display	Controller failure	Y
Ambient temperature	10~32°C	External material	SUS 304	Door ajar	Y
Noise (db)	55	Compressor brand/QTY	Secop/1	Power failure backup	8 hours
Electrical Data		Dimensions		Other	
Voltage frequency( V/Hz)	220V 50Hz/60Hz	Package size(W*D*H)(mm)	1880*820*2150	Inside door	12 shelves
Power (W)	1040W	Capacity (L/cu.ft)	1500/53	Remote alarm port	Optional
Current (A)	5.4A	Exterior size(W*D*H)(mm)	1800*775*1965	Access port QTY/Dia.	Optional
Power consumption	8kWh/24h	Interior size(W*D*H)(mm)	1680*595*1312	Internal material	SUS 304