

Large Capacity Refrigerated centrifuge DL-6MC

Feature

1. With the self-diagnostic and protection system for imbalance, over-speed, over-temperature. Greatly ensured the safe of users and machine.
2. Adopt microprocessor control, high efficiency frequency inverter motor, R404a environmentally friendly refrigerant compressor.
3. Emergency door lid release system, in case of unexpected power off.
4. Equipped with wheels, easy to short distance move.
5. Both touch and button operation, high-light touch LCD display shows full parameters.
6. 35 programs can be editable, storable and retrievable, users can set 9 programs for daily use,




Parameters

Model	DL-6MC	Max Speed	6000r/min
Max RCF	6680xg	Max Capacity	6×1000ml
Speed Accuracy	±10r/min	Timer Range	0-99h 59 min
Temperature Range	-20~+40°C	control&drive	microprocessor control/direct drive
Acc/dec Rate	0-9(Grade)	Program Store	35
Motor Power	3.8KW	Refrigerator Power	2.2KW
Noise	≤60dB(A)	Power Supply	AC220V /380V 50-60Hz 35A
Weight(Without Rotor)	260kg	Dimensions	730×840×1200mm(L×W×H)
Refrigeration System	Big Brand High-performance Compressor CFC-free Refrigeration System		
Rotor	No.1 AR 6*500ml (6000rpm 6680xg) No.2 SR 6*1000ml Round Cup (4200rom 5600xg) No.3 SR 6*1000ml Oval Cup (4000rom 5600xg)		

Application:



Blood Bank Center

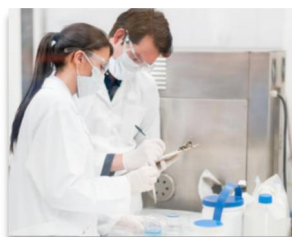

Bio-Technology Company



Chemistry Lab



Bio-Pharmaceutical